

Product specifications

Name	Anti-h NT-proBNP RC1306 SPTN-5
Specificity	Antibody recognizes human N-terminal pro-B-type natriuretic peptide
Description	Recombinant chimeric antibody containing variable region from mouse IgG and constant region from human IgG ₁ , cultured <i>in vitro</i> under conditions free from animal-derived components.
Product code	140010
Product buffer solution	50 mM Na-citrate, pH 6.0, 0.9 % NaCl, 0.095 % NaN ₃ as a preservative
Shelf life and storage	Unspecified, storage at 2–8 °C
Subclass	human IgG ₁
Analyte description	The plasma concentrations of both the N-terminal fragment of Pro-B-type natriuretic peptide (NT-proBNP) and prohormone of brain natriuretic peptide (proBNP) are significantly increased in patients with asymptomatic and symptomatic left ventricular dysfunction.

Parameters tested on each lot


Product appearance	Liquid, may turn slightly opaque during storage
Product concentration	5.0 mg/ml (+/- 10 %) (A280 nm, 1 mg/ml, 1 cm = 1.6)
Immunoreactivity	80–120 % compared to the reference sample in an FIA test
IEF Profile	8.1–8.9
Purity	≥ 95 %

Kinetic parameters

Association rate constant	Not Determined (N/D)
Dissociation rate constant	N/D
Affinity constant	N/D
Determination method	-
Determination antigen	-

Legal disclaimer

Medix Biochemica products meet their specifications if transported, stored and used according to the instructions. MedixMAB and MedixAntigens are trademarks of Medix Biochemica and may not be used or reproduced without Medix Biochemica's written permission.

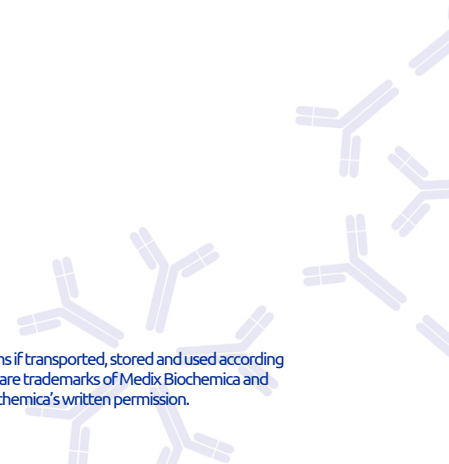

Cross-reactivities Does not recognize canine NT-proBNP (HyTest, 8CNT9).

Epitope Amino acids 12-28 (LETSGLQEQRNHLQGKL)

Pair recommendations

		DETECTION						
		1306/ RC1306	1307/ RC1307	1308	1309	1310	1311	1312
CAPTURE	1306 / RC1306	-	+	-	+	+	+	+
	1307/ RC1307	+	-	+	+	+	+	+
	1308	-	+	-	+	+	+	+
	1309	+	+	+	-	-	-	-
	1310	+	+	+	-	-	-	-
	1311	+	+	+	-	-	-	-
	1312	+	+	+	-	-	-	-

Please note that pair recommendations are based on results obtained by our laboratory. Equally good results may be obtained using other pairs and therefore these recommendations are only indicative.

Platforms tested FIA

Antigens tested Recombinant NT-proBNP antigen, Medix Biochemica 610090 and recombinant proBNP antigen Medix Biochemica 710017, 710042 and 710043.

TEMPERATURE, TIME	RESULT
-70 °C, 21 days	OK
-20 °C, 21 days	OK
+4 °C, 21 days	OK
+35 °C, 21 days	OK
+45 °C, 7 days	OK

Stability testing is performed in the product buffer to see whether different temperatures affect the antigen binding, charge or composition of the antibody. Please note that the shelf life given on the first page is based on real time stability testing at 2–8 °C in the product buffer.

Miscellaneous -

References -

